

SC\_BODY\_FAT\_B


1:

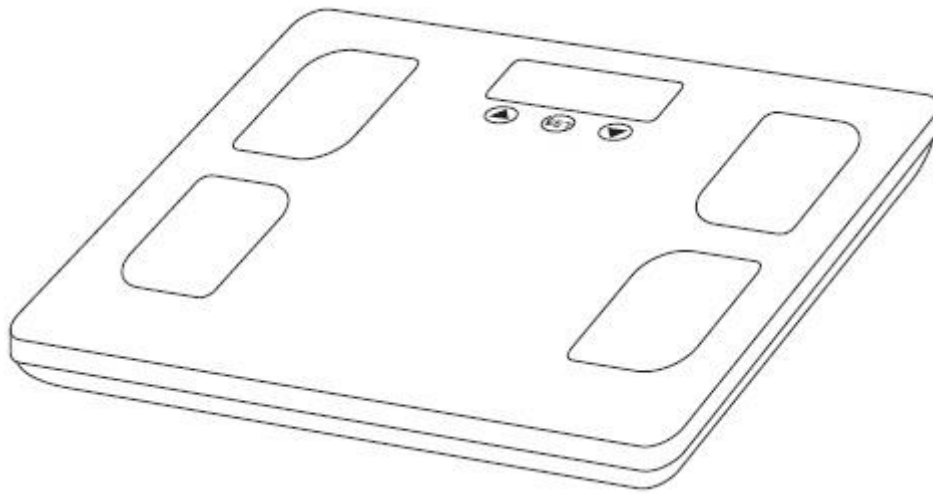
2:

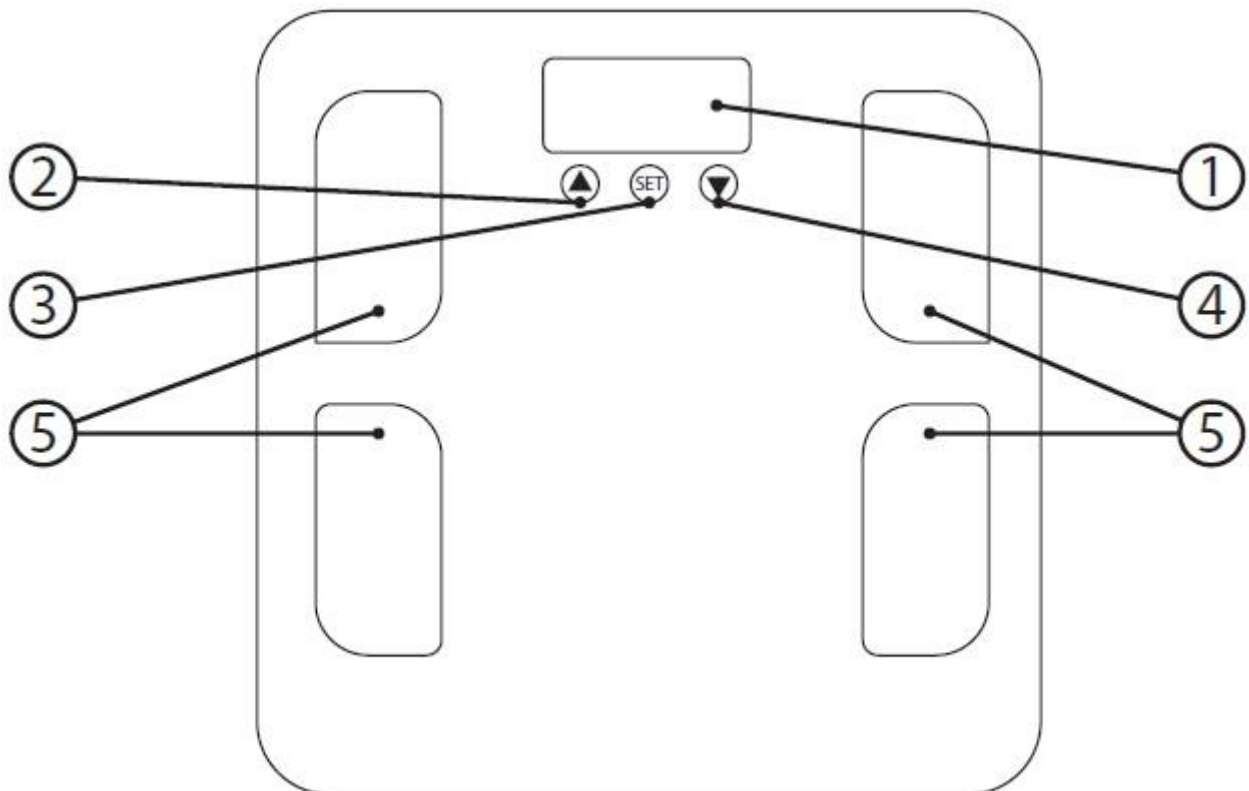
3:

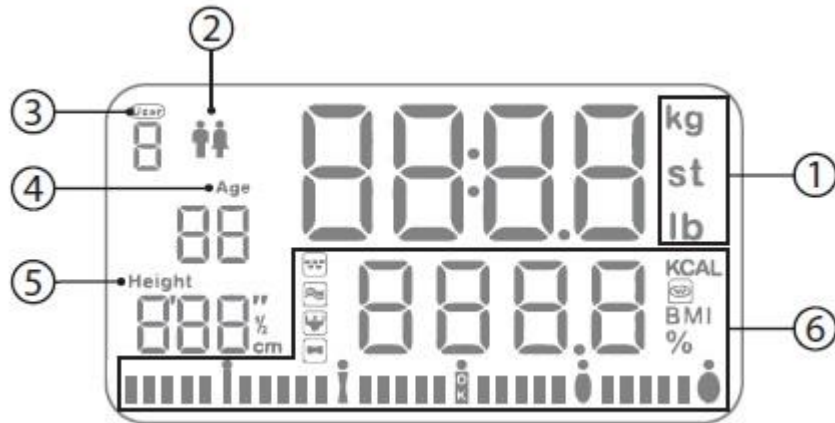4:

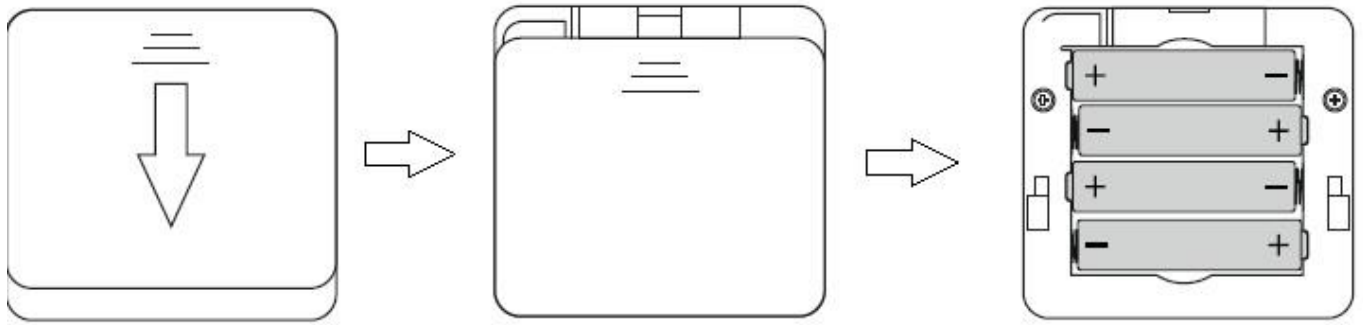5:

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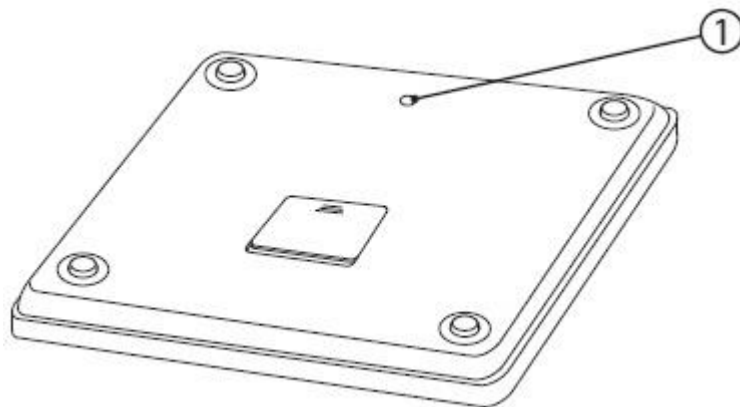WEEE-Reg.-Nr.: DE 78960174

### Manual

Thank you for choosing a product of SEBSON. Read the instructions carefully and keep them for looking up before start-up this product. A long and failure-free operation can only be guaranteed if the installation and start-up is done properly. Before starting up this product, please check for external damage. In case of visible external damage we immediately dis-advise the start-up of this product! Please remove all packaging material of this product before the first start-up and check the content for completeness!

### Intended use:

This body fat scale of the brand SEBSON you can determine besides the body weight further body values like body fat, muscle mass etc.

**Scope of delivery:**

- SC\_BODY\_FAT\_B

**General safety information:**

**Danger of suffocation:** Keep this product away from children and animals! The packaging contains small parts that can be swallowed! Keep the device away from children under the age of 8!

**Note:** This device is not a medical device and therefore does not replace a medical measure!

- This device is for private use only. This device is not suitable for commercial or medical use! This device is not officially calibrated!
- Persons with reduced physical, sensory or mental abilities or lack of knowledge and experience may only put this device into operation if they are able to handle it safely and are aware of the dangers in handling this device!
- The remodeling of this product is not allowed! The proper and safe functioning of the device can no longer be guaranteed if it is modified! In addition, this will invalidate the warranty claim!
- This product is only suited for indoor-use and not for outdoor-use!
- Do not use the device near high-frequency electromagnetic transmitters. The transmitters may cause interference to this equipment!
- Protect the device from shocks and vibrations.
- Avoid overloading the device! Observe the maximum load capacity specified in the "Technical data". Overloading can damage the device.
- Please do not step on the scale with wet feet. There is a danger of slipping. The scale may also be damaged by moisture!
- Please do not store the device vertically!
- Do not drop the device because this can damage it! Don't drop anything on the device! **Note:** There is a danger of glass splinters here!
- The manufacture does not accept any responsibility for injuries or damage resulting from improper use of this product!
- Any other use than the description may result in damage to the product. So the dangers are such as electric shock, short circuit, fire etc.
- We reserve the right to make changes to the product, the packaging or the operating instructions as well as printing errors.
- Protect the product from external damage, heavy dirt, moisture and excessive heat!
- In order that the warranty claim does not expire please use only the supplied items. Repairs may only be done by an expert and only by using original spare parts!

**General information about the product:**

Following values can be determined with this scale:

- Weight (kg/ lb/ st)
- BMI
- Body fat rate (%)
- Body water (%)
- Muscle mass (%)
- Bone mass (%)
- BMR (kcal)
- Visceral fat

Up to 8 profiles of different people can be created and analyzed.

**Product overview (Fig.2):**

1. Display
2. Set values (increase)
3. Switch on and enter key
4. Set values (decrease)
5. Measurement sensors

**Display indication (Fig.3):**

1. Weight units (kg/ st/ lb)
2. Gender
3. User
4. Age
5. Body size
6. Determined body measurements

**Explanation of symbols:**

	Body fat
	Body hydration
	Body muscle
	Body bone
	Visceral fat
	BMI-Bar graph
<b>BMI</b>	Body mass index (BMI)
<b>KCAL</b>	Calorie requirement

**Groups of persons who are not suitable for the body analysis function:**

To analyze the body, the balance sends a weak electrical impulse through the body. The following groups of people should not use this device for body analysis:

- Children under the age of 10
- Persons older than 70 years of age
- pregnant women
- People with pacemakers or other medical implants

This can lead to incorrect measured values for the following groups of persons:

- dialysis patients
- Persons with high fever
- Persons with severely swollen legs or oedemas
- Persons whose bodies have high water retention or heavily drained persons
- Persons suffering from osteoporosis
- Competitive athletes and bodybuilders
- Teenagers under 17 years of age
- Persons who spend more than 10 hours per week intensively doing sports

- Persons with a lot of cornea under their feet

**Inserting or changing the batteries (Fig.4):**

- **Note:** Please observe the general instructions for using batteries!
- Insert 4 LR03 AAA/1.5V batteries into the battery compartment. **Note:** Make sure that the polarity is correct when inserting the batteries!
- Close the battery compartment again after inserting the battery. Please make sure that the battery compartment is always closed! Do not use the device if the cover of the battery compartment is damaged or missing!
- Replace the batteries as soon as the symbol "Lo" appears in the display or the device cannot be switched on at all!

**General instructions for using batteries:**

**Danger:** Never expose batteries to extreme conditions, such as heat or strong sunlight! There is an acute risk of explosion and leakage!

**Danger:** Make sure that batteries can not be swallowed by children, animals and adults! In particular, keep children and animals away from this product! If one or more batteries are swallowed, seek medical attention immediately! Here is an acute danger to life. The substances contained in the battery can cause severe poisoning!

**Caution:** Remove leaking batteries from the device immediately! The battery acid contained in the battery can damage the unit! Avoid direct skin contact with leaking batteries! Wash your hands immediately after contact and take care that the battery acid does not get into your skin, eyes or mucous membranes! Dispose of leaking and old batteries according to the regulations.

- Replace old batteries with batteries of the same type!
- Do not use rechargeable batteries. They have a too low voltage!
- Remove the batteries from the battery compartment when you are not using the remote control for more than a month. Make sure that the battery compartment is always closed when the battery is inserted!
- Only operate the device when the battery compartment is closed!
- Non-rechargeable batteries are not to be re-charged!
- New batteries should not be used together with used batteries!
- Always use batteries with the same type!
- Never short the terminals of the batteries!

**Change unit (Fig.5):**

**1:** The button on the back of the scale can be used to switch on or change the weight unit by lightly touching it with your finger.

**Initial start up:**

A short warm-up phase is recommended for initial commissioning or when relocating.

Please follow the procedure as described:

1. Place the device on a firm, vibration-free and level surface. **Note:** Please make sure that all 4 leveling feet of the balance touch the floor.
2. Place one foot on the scale to activate the device.
3. Remove your foot and wait for the scale to turn off. Leave the scale turned off for about 10 seconds.
4. The scale is now ready for operation.

### Measure weight:

1. Place the device on a firm, vibration-free and level surface. Please make sure that all 4 leveling feet of the balance touch the floor.
2. If possible, stand in the middle of the scale. Do not stand on the display! **Attention:** If you stand on the edge of the scale, the scale may tilt to the side. There is a risk of injury!
3. The scale is now activated and measures your weight.
4. During the measurement, hold as still as possible to avoid measuring errors.
5. The measurement is finished as soon as the weight display starts flashing.
6. *Automatic switch-off function:* After approx. 8 seconds of inactivity, the device switches itself off automatically.

### Store personal data for body analysis:

In order to determine, for example, the percentage of body fat or body water of a person with this device, information such as age, sex and height of this person must first be entered into the device. Up to 8 persons can be stored. To do this, please proceed as described:

1. Press SET-key
2. Select "User" with the ▲ and ▼ keys and confirm with "SET".
3. Select gender  with the ▲ and ▼ keys and confirm with "SET".
4. Select "Age" with the ▲ and ▼ keys and confirm with "SET".
5. Select "Height" with the ▲ and ▼ keys and confirm with "SET".

### Perform body analysis:

It is recommended to always perform the body analysis at the same time of day with the water balance remaining as constant as possible (e.g. in the morning after the first visit to the toilet). Wear no or only very light clothing during the measurement. **Note:** The body analysis can only be carried out barefoot and only standing up! Please enter the scale with clean and dry feet!

The values are determined as described:

1. Place the device on a firm, vibration-free and level surface. Please make sure that all 4 leveling feet of the balance touch the floor.
2. Press SET-key
3. Select "User" with the ▲ and ▼ keys and confirm with "SET".
4. Confirm gender  with "SET".
5. Confirm "Age" with "SET".
6. Confirm "Height" with "SET".
7. As soon as "0.0" appears in the display, the measurement can be started.
8. Place yourself on the metal surfaces of the balance with your legs as straight as possible. Your knees must not touch during the measurement. Stand as still as possible during the measurement.
9. Your weight is shown on the display first. Stay still on the scale. The weight display will stop flashing after a few seconds. Then the other values such as the percentage of body fat etc. are determined. This is indicated on the display by the following symbol: . Continue to stand still on the scale.
10. After a successful measurement, the following values are displayed one after the other below the weight display in the following sequence: BMI → body fat content → body water content → muscle mass content → bone mass content → calorie requirement → visceral fat content
11. Then exit the balance. The measured values are displayed three times in succession. The balance then switches itself off automatically. The values are saved.

**Display stored values:**

1. Press SET-key
2. Select "User" with the ▲ and ▼ keys and confirm with "SET".
3. Confirm gender  with "SET".
4. Confirm "Age" with "SET".
5. Confirm "Height" with "SET".
6. As soon as "0.0" appears on the display, press this key: ▲. The values determined last are then shown three times in succession in the display. The unit then switches itself off automatically.

**Delete stored values:**

1. Press SET-key
2. Select the "User" to be deleted with the ▲ and ▼ keys and confirm with "SET".
3. Confirm gender  with "SET".
4. Confirm "Age" with "SET".
5. Select "Height" with ▲ >250 or with ▼ <100 until 0 is displayed and confirm with "SET". The stored data for this user have been deleted.

**Evaluation table:**

Age						KCAL	
<20	low	<12		<46		<1.200	
	normal	12~20	46~69	38~56	2.0~4.1	1.200~2.600	1~9
	high	>20	>69	>56		>2.600	<10
20-30	low	<13	<48	<43		<1.150	
	normal	13~21	48~69	43~64	2.1~4.0	1.150~2.800	1~9
	high	>21	>69	>64		>2.800	>10
30-40	low	<15	<45	<48		<1.100	
	normal	15~23	45~68	48~69	1.8~4.0	1.100~2.400	1~9
	high	>23	>68	>69		<2.400	>10
40-50	low	<16	<43	<46		<1.030	
	normal	16~24	43~68	46~62	1.9~3.8	1.030~2.150	1~9
	high	>24	>68	>62		<2.150	>10
50-60	low	<17	<43	<44		<980	
	normal	17~25	43~65	44~59	1.9~3.7	980~2.060	1~9
	high	>25	>65	<59		<2.060	>10
>60	low	<18	<42	<41		<720	
	normal	18~26	42~65	41~58	1.6~3.7	720~1.750	1~9
	high	>26	>65	<58		<1.750	>10

**Table for the Body Mass Index:**

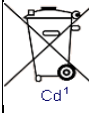
					
Gender:	Underweight	Normal weight	Overweight	Obesity	Severe obesity
	<19	19-24	25-30	31-40	>40

	<20	20-25	26-30	31-40	>40
---	-----	-------	-------	-------	-----

Solve problems:	
	Replace the batteries!
	Max. load capacity has been exceeded. Please leave the scale immediately to avoid damaging the device!
	If this error message appears on the display during the measuring process, exit the balance and restart the device.

Technical Data:	
Operating voltage:	6V DC
Battery type:	4 x 1,5V AAA
Measuring units:	kg/ lb/ st
Minimum load:	≥5kg/ 11lb
Maximum resilience:	≤180kg/ 400lb
Steps:	100g/ 0,2lb

	<b>Disposal note:</b>
<p>This product may not be disposed of with household waste. At the end of its useful life, the product must be returned to a point of acceptance for the recycling of electric and electronic device. Please ask your local municipality for the point of acceptance! As the devices contain valuable recyclable materials, they should be sent for recycling, in order not to harm the environment or human health from uncontrolled waste disposal.</p>	

	<b>Proper disposal of batteries and accumulators:</b>
<p>Batteries and rechargeable batteries may not be disposed in the household waste. The end-consumer are legally obligated to the return of all used batteries and/or accumulators. After usage you can return batteries and / or accumulators gratuitous in the sales point or in their direct vicinity (for example at municipal collection points and trade). Otherwise, you can also return the batteries and / or accumulators by post to us.</p>	
  	<p>Batteries and / or accumulators which contain harmful substances are marked with the symbol of a crossed-out waste bin and may not be disposed of with household waste. Under the dustbin symbol is the chemical name of the harmful substance - in the example below - "CD" for cadmium, "Pb" stands for lead, "HG" for mercury.</p>

If not disposed of properly, the substances contained in the batteries can cause health and environmental hazards. These dangers are averted by the separate collection and proper disposal of used batteries.

#### Cleaning:

Use a slightly moistened cloth for cleaning. In order not to damage the product, it is essential to avoid water penetrating into the housing! Never clean the device under running water! Do not use solvent-containing cleaner for cleaning. Do not use solvent-containing cleaner for cleaning!

#### CE Declaration of Conformity:

This product has been prepared according to the listed below guidelines:

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